


GAS SELECTION PANEL

The Gas Selection Panel, when used in conjunction with the Sechrist Monoplace Hyperbaric Chamber system, will allow the chamber operator to select either oxygen or medical grade air to be supplied to the chamber. The panel incorporates pressure gauges for the oxygen, medical grade air and the selected gas source. The design of the gas selector valve utilizes a “shear type” sealing method, and the valve vents to ambient.

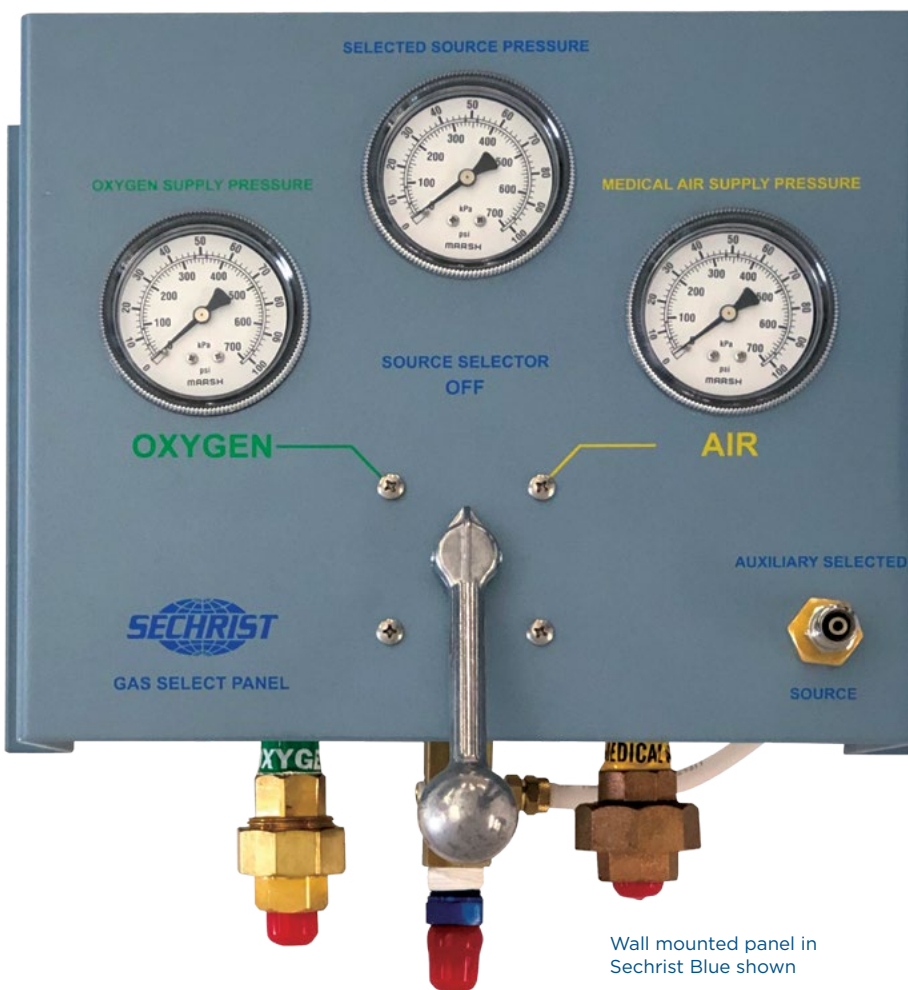

Wall mounted panel in Sechrist Blue shown

PRODUCT	PART NUMBER
Gas Selection Panel	
Sechrist Blue	20200-00
Diamond White	20200-01
Sahara Beige	20200-02
Platinum	20200-03
Crystal Green	20200-04
Jade Green	20200-05
Pressure Gauge	27010
Selector Valve	50023
Supply Hose Adapter	HB 701-02

PRODUCT	PART NUMBER
Height	13.5 inches (34.29 cm)
Width	14 inches (35.56 cm)
Depth	7 inches (17.78 cm)
Operating Pressure	50 to 70 p.s.i. (3.52 to 4.92 kg/cm ²)

FOR COMPLETE DOCUMENTATION, SEE PACKAGE LABELING ACCOMPANYING DEVICE